




LED Silicone Beauty Mask

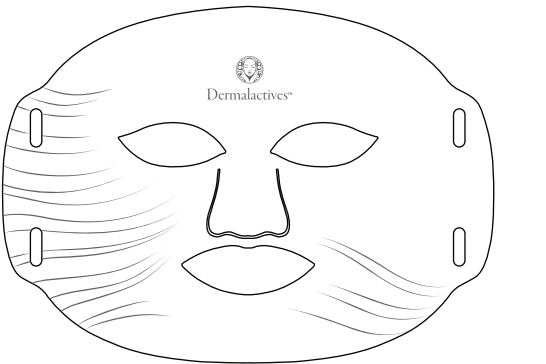


User manual

Please read the manual carefully before using the product
Please keep it in case need

Experience professional-inspired LED light therapy designed to complement your skincare routine. The Dermalactives LED Silicone Beauty Mask combines flexible silicone comfort with multiple light wavelengths to help promote healthier-looking, more radiant skin—all from the comfort of your home.

►► INCLUDING

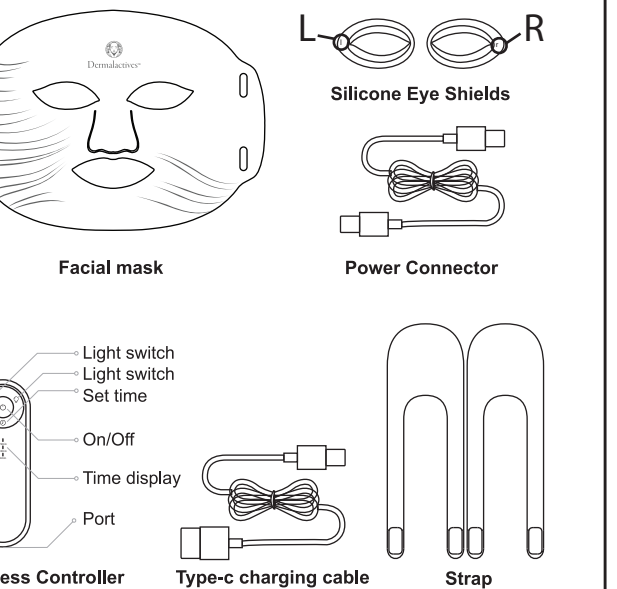

►► SPECIFICATION

Model	D-0057S
Battery	800 mAh Lithium-Ion
Charging Port	USB Type-C
Charging Voltage	DC 5V
Power	3W
Operating Voltage	3.7V
Weight	162 g

LED LIGHT GUIDE

Color	Benefit
Red	Helps support firmer-looking skin and promotes a more youthful appearance
Blue	Helps improve the appearance of blemish-prone skin
Green	Helps promote a more even-looking skin tone
Yellow	Helps improve the appearance of dull-looking skin
Purple	Combination treatment
Cyan	Helps soothe the appearance of stressed skin
White	Multi-spectrum cosmetic treatment

SUGGESTED WEEKLY ROUTINE

Goal	Frequency
General Skin Maintenance	3 times per week
Dull-Looking Skin	4–5 times per week
Blemish-Prone Skin	4–5 times per week
Mature Skin	3–5 times per week

►► HOW TO USE

Step 1	Prepare the Mask: Attach the silicone eye shields before use. Position the shield marked "L" over the left eye opening and the shield marked "R" over the right eye opening, ensuring they are securely aligned.
Step 2	Wear the Mask: Attach and adjust the head straps for a comfortable, secure fit. Place the mask over your face, then press and hold the Power button for approximately 1.5 seconds to turn the device on.
Step 3	Select Your Treatment: Press the Light Mode button to choose your desired LED light color. Then press the Timer button to select your treatment duration: • 5 minutes • 10 minutes • 15 minutes The device will begin treatment automatically.
Step 4	Complete Your Routine: Once the selected treatment time has ended, the device will automatically power off. Remove the mask and continue with your preferred skincare routine.
Remark	Recommended Treatment Schedule For best cosmetic results, use the mask 3–5 times per week for 4 consecutive weeks. After completing one cycle, continue using the device as needed to maintain results.

Before Your First Treatment

Perform a skin sensitivity test before using the device for the first time.

1. Place the mask against the inside of your forearm.
2. Turn the device on and select a 10-minute treatment.
3. Wait 24 hours.
4. If redness, itching, or irritation occurs, discontinue use and consult a healthcare professional before using the device again.

Before Each Use

- Thoroughly cleanse and dry your skin.
- Remove all makeup, sunscreen, and skincare products before treatment.
- Ensure the mask and eye shields are clean before each use.

After Each Use

Mild redness immediately after treatment is normal and typically subsides within **1–2 hours**. If you experience excessive warmth, burning, pain, or prolonged irritation, discontinue use and consult a healthcare professional.


►► CONTROLLER FUNCTIONS

Power

Press and hold the **Power** button for approximately **1.5 seconds** to turn the device on or off.
Default setting: Red LED mode.

Light

While the device is on, press the **Light Mode** button to cycle through the available LED light colors.
If your controller includes two navigation buttons:
Right button: Scroll forward through the light modes.
Left button: Scroll backward through the light modes.

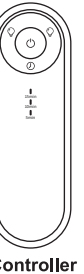

Timer

Press the **Timer** button to select your treatment duration:

- **5 minutes** (default)
- **10 minutes**
- **15 minutes**

The device will automatically shut off when the selected treatment is complete.

►► SAFETY PRECAUTIONS

Please read the following instructions carefully before using this device.

1. For indoor use only. This device is not waterproof. Do not immerse it in water or expose it to excessive moisture. Clean only as instructed.
2. Store the device in a cool, dry place away from direct sunlight, excessive heat, and humidity.
3. If the device is not used for an extended period, recharge the battery at least once every three months to help maintain battery performance.
4. Do not disassemble, modify, or attempt to repair the device yourself. Doing so may damage the product and void the warranty.
5. Use this device only as directed in this manual. Do not exceed the recommended treatment frequency or duration.
6. If all three battery indicator lights flash simultaneously when the controller is turned on, the battery is low and should be recharged before use.
7. This device is not intended for use by children or by individuals with reduced physical, sensory, or mental capabilities unless they are supervised by a responsible adult.
8. Discontinue use immediately if you experience excessive discomfort, prolonged redness, burning, swelling, or any unexpected skin reaction. If symptoms persist, consult a healthcare professional.

9. Avoid looking directly at the LED lights during operation. Always install the supplied eye shields before facial treatments.
10. Use only the supplied charging cable and approved accessories.

Do Not Use This Device If You:

- Are pregnant, planning to become pregnant, or breastfeeding.
- Have a known sensitivity to light or are taking medications that may increase photosensitivity.
- Have epilepsy or a history of seizures triggered by flashing lights.
- Have active skin cancer or suspicious skin lesions in the treatment area.
- Have lupus, albinism, photosensitivity eczema, or other photosensitive medical conditions.
- Have open wounds, severe acne lesions, active skin infections, or significant irritation in the treatment area.
- Have undergone recent facial surgery, laser treatments, chemical peels, or other invasive cosmetic procedures unless approved by your physician.
- Have implanted electronic medical devices, such as a pacemaker or defibrillator, unless cleared by your physician. Important: This device is intended for cosmetic use only. It is not intended to diagnose, treat, cure, or prevent any disease. If you have a medical condition or concerns about using LED light therapy, consult your physician before use.

Charging Instructions

Before using the device for the first time, fully charge the controller.
1. Connect the supplied **USB Type-C charging cable** to the controller.
2. Connect the cable to a certified USB power adapter (5V DC).
3. The battery indicator will **flash while charging**.
4. When the indicator remains **solid**, the battery is fully charged.
5. Disconnect the charging cable before using the device.

Important

- Do not operate the device while it is charging.
- Use only the supplied USB Type-C cable or a compatible replacement.
- Do not charge the device near water or in excessively humid environments.

►► CARE & MAINTENANCE

Proper care will help maximize the performance and lifespan of your device.

Cleaning

- Turn the device off before cleaning.
- Disconnect the controller from the mask.
- Wipe the silicone surface with a soft, slightly damp cloth.
- Dry thoroughly before storing.
- Do not immerse the mask, controller, or charging cable in water.
- Do not use alcohol, acetone, abrasive cleaners, or harsh chemicals.

Storage

- Store in a cool, dry place at room temperature.
- Avoid direct sunlight, excessive heat, humidity, and dusty environments.
- Do not place heavy objects on the mask.
- Recharge the battery at least once every three months if the device is not in regular use.

Battery Disposal

This device contains a rechargeable **lithium-ion battery**.

Do not dispose of the device with household waste. Follow your local regulations for recycling electronic devices and rechargeable batteries.

►► LIMITED WARRANTY

Dermalactives warrants the controller against defects in materials and workmanship for **one (1) year** from the original date of purchase when purchased through an authorized retailer.

This warranty does not cover:

- Normal wear and cosmetic damage
- Damage caused by misuse, abuse, neglect, or improper storage
- Damage resulting from unauthorized repair, modification, or disassembly
- Damage caused by accidents, impact, or dropping the device
- Consumable accessories, including straps, eye shields, charging cable, or other replaceable components

For warranty assistance, replacement parts, or product support, please contact Dermalactives Customer Support or visit www.dermalactives.com.

* The right to interpret this warranty card belongs to the brand dealer or manufacturer

Product name			
Product number		Production Date	
Customer Name			
Phone number		Customer Fax	
Problem			

