



Cautions

- Do Not Use this product if you have broken skin, subcutaneous implants, undergone facial plastic surgery, or are pregnant or menstruating.
- It is only intended for facial and neck care. Do Not Use on other parts of the body.
- Avoid contact with pimples, acne, or rosacea.
- People with sensitive skin are not recommended to use it daily; once a week is best.


Charging & Maintenance

- The main unit has a built-in rechargeable battery. Charge the device using the included USB Type-C charging cable.
- The device cannot be powered on while charging. Disconnect the charger once the device is fully charged to help preserve battery life.
- If the device will not be used for an extended period, recharge it periodically to maintain battery performance.
- Disconnect the device from the charger before cleaning or performing maintenance. Do not immerse or wash the entire unit in water. For daily cleaning, wipe the device with a lint-free wipe or a soft cloth lightly dampened with a mild cleaning solution.
- Do not store or use the device in humid environments, such as bathrooms, as moisture may cause electrical damage or a short circuit.


Parameters

Product Name	V-Lift Sculpt
Material	ABS
Rated Power	2.5W
Battery	3.7 V/800mAh
Voltage	5V—1A
Product Size	128*62*40mm
Weight	110.9g
Charging Time	2-3 hours

Accessories

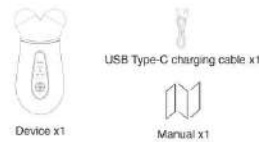

Dermalactives™
ELASTIN & COLLAGEN REJUVENATION

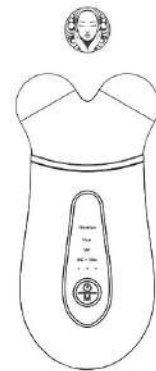

V-Lift Sculpt
User Manual

Functions

Thank you for choosing the V-Lift Sculpt which combines vibration, heat therapy, and microcurrent technology. To ensure safe and effective use of this product, please read this manual carefully before use and keep it for future reference.

Vibration Mode

High-Frequency Vibration relaxes facial and neck muscles, shapes facial contours, and boosts the absorption of skincare products.

Heat Mode

Heat Therapy opens pores, soothes skin, and promotes blood circulation.


Microcurrent Mode (MC)

Microcurrent stimulation activates muscles, tightens saggy skin, and promotes a firmer appearance.

Microcurrent + Vibration Mode (MC + Vibe)

Improves lymphatic drainage, reduces edema, and helps draw out impurities.

Product Structure

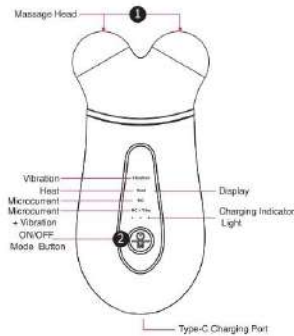

1 Massage Head

Directly contacts the skin to deliver warmth, vibration, and microcurrent energy.

2 ON/OFF & Mode Button

Press and hold for 1.5 seconds to power on/off; short press to switch working modes.

How to Use


Preparation Before Use

- Skin Preparation: Ensure your face is clean and free of oil, dirt, and makeup before use.
- When using the microcurrent mode, apply a water-based gel, essence, serum, or facial mask. Do Not Use on dry skin.
- Device Preparation: Fully charge the device before using it for the first time. The indicator light will illuminate while charging and remain solid when the device is fully charged.
- Safety Check: Inspect the whole device for cracks or abnormalities and if the buttons are responsive.

Operation Instructions

- Press and hold the power button for 1.5 seconds to turn on the device. The device will automatically start in Vibration Mode.
- Each short press switches to the next mode, with the corresponding mode icon lighting up on the display area.
- Press and hold the button for 1.5 seconds to power off.
- Recommended Usage Duration: Limit each full facial treatment to 10–15 minutes. Use the device 3–5 times per week to avoid over-treating your skin.

Treatment Area

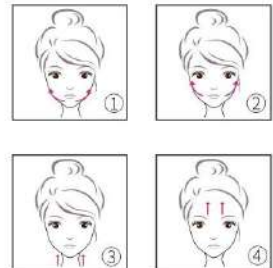

- ① Lift upward from the chin to the cheeks.
- ② Lift from the cheeks toward the outer corners of the eyes.
- ③ Lift upward from the lower neck to the chin.
- ④ Lift upward from the lower forehead to the upper forehead.